

GEMINI

COMPACT MICROPLATE
PROCESSOR



TUBE IN – RESULT OUT

FULLY AUTOMATED MICROPLATE PROCESSOR FOR LOW THROUGHPUT APPLICATIONS

- Total walk away system
- Flexible time management with LIS connectivity
- Modularity – open system
- Ease of use, installation and maintenance
- Assay performance with disposable tips
- High process safety with triple pipetting control
- Convenient programming of assay protocols
- Proven reliability
- The system is in compliance with the In-vitro diagnostic directive 98/79/EC

Reading

- Absorbance reader with up to 8 different filters
- Single or dual wavelength reading
- End point and kinetics
- Shaking

Incubation


- Two independently controlled incubator chambers (up to 45 °C) with shaking (20 Hz)
- Three room temperature incubation chambers

Washing

- 2 x 8-channel washer (aspirate and dispense)
- Three wash buffers
- Sweep mode, soak, top and bottom wash, variable pump speed

Capacity and throughput

- Initial loading of up to 192 samples
- Continuous loading
- Two deep-well plate positions for predilution and/or archiving


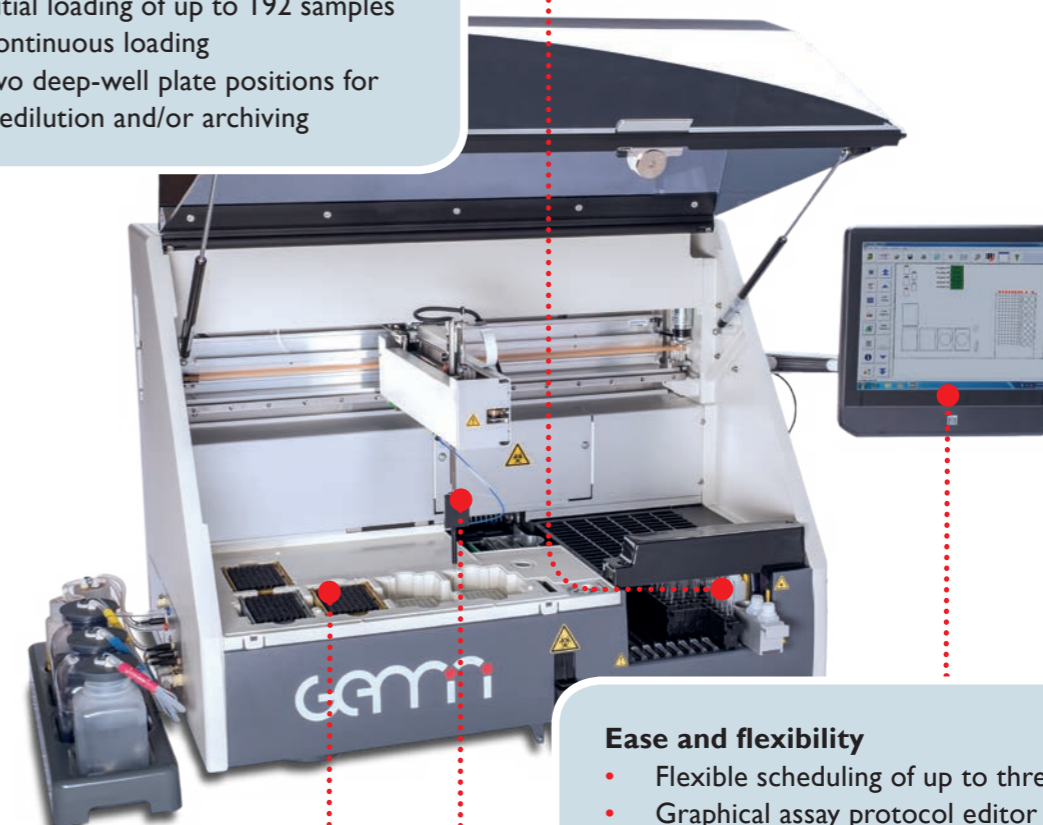


Ease and flexibility

- Flexible scheduling of up to three plates
- Graphical assay protocol editor
- Combination of multiple assays on one plate
- Proven software reliability
- Bi-directional LIS connectivity, ideal with random access instruments
- Based on established EIA software platform with more than 4,500 placements worldwide

Safety and in-process control

- Pipetting pressure monitoring
- Capacitive liquid level detection for optimal identification of fibrin clots and insufficient volume
- Colorimetric dispense verification
- Positive identification of barcoded samples and reagents

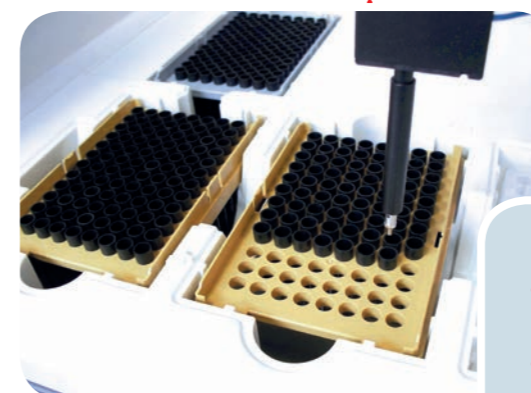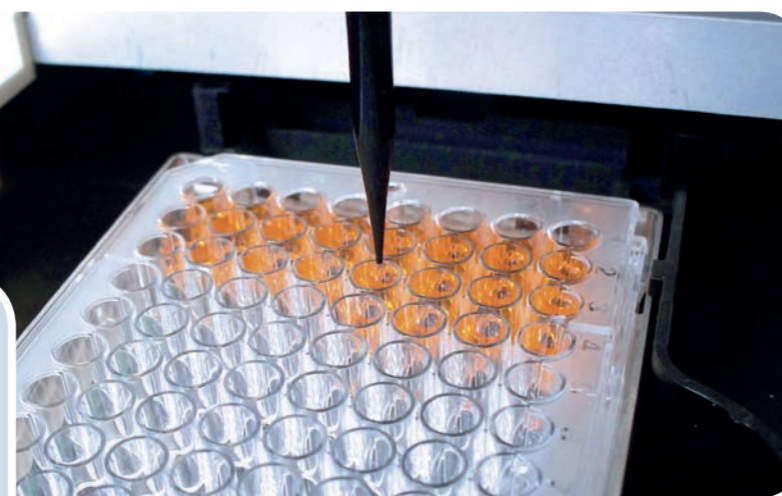

Pipetting

- Liquid pipettor for disposable tips (300 µl or 1100 µl)
- Flexibly programmable for single pipetting, multiple pipetting, multiple dilution steps, mixing and archiving

SPECIFICATIONS

READER

- Photometric range: 0 to 3.0 OD
- Spectral range: 400 nm to 700 nm (up to 8 filters)
- Read time: < 15 sec single, < 30 sec dual
- Precision: 1 % CV at 1.0 OD
- Accuracy: ± 0.005 OD or 2.5 % (whichever is greater)
- Linearity: 0 to 2.000 OD ± 1 %

PIPETTOR

- Min / max. volumes: 10 μ l to 300 μ l with 300 μ l tip
301 μ l to 1000 μ l with 1100 μ l tip
- Precision (single dispense): < 3 % CV at 20 μ l
< 3 % CV at 100 μ l
- Precision (multi dispense): < 10 % CV at 16 x 20 μ l
< 3 % CV at 8 x 100 μ l
- Features: Pipetting pressure monitoring, capacitive liquid level detection, tip detection, mixing, multiple dilution steps, archiving

CAPACITY

- Sample and reagent capacity: Up to 192 samples
Flexible: e.g. 144 samples + 8 reagents + 16 controls

INCUBATOR

- Temperature range: Up to 45 °C
- Temperature uniformity: ± 1.5 °C (with in-process temperature monitoring)
- Shaking: 20 Hz

All values are achieved under optimal conditions and can vary depending on environmental conditions, instrument status and processing conditions.

WASHER

- Precision: 10 % CV at 300 μ l
- Residual volume: < 2.5 μ l in U-bottom (mean)
< 4 μ l in flat bottom (mean)
- Wash buffer capacity: 3 wash buffers
- Modes: Sweep mode, soak, top and bottom wash, variable pump speed

DIMENSIONS

- Width x Depth x Height
- 125 cm x 65 cm* x 75 cm**
- 49" x 26"* x 30"***

EASY INTEGRATION – NO LIMITS

Additional functionality through middleware software

- Easy integration: Wide range of interfaces allows consolidation of results from other instruments
- Connects host: Provides smooth real time bi-directional communication between device and LIS
- Connects instruments: Operates with single or multiple GEMINI installations

* w/o Waste Bag

** Cover Closed

diatron●●

We are a member of the STRATEC Group | www.diatron.com | sales@diatron.com

Diatron MI Zrt makes best efforts to provide accurate information and specifications herein, whereby it may be possible for omissions or errors to have occurred. Diatron MI Zrt shall have no liability for the accuracy and completeness of the information and cannot be held liable for any representations or warranties, expressed or implied, and any third party claims or losses of any damages. Diatron MI Zrt reserves the right to make any changes to the information, including specifications, without notice. The applications and products herein may not be available in all markets. ©2019, Diatron MI Zrt, Hungary, all rights reserved.

Version Nr. GEMINI I 90605V2

